



® Knowledge Beyond Measure.

# Portable Test Aerosol Generator

Models 3073



## Provides stable test aerosol for demanding applications

The Portable Test Aerosol Generator 3073 is a truly portable high-end test aerosol generator. Polydisperse test aerosols are generated by atomizing various oils or salt solutions where the mode of the produced particle size distribution is typically between 0.15 and 0.3  $\mu\text{m}$ . Monodisperse particles such as PSL can be generated as well from liquid suspensions. The innovative new flow control via nozzle pressure sensor and proportional valve achieves highly stable low particle production rates. Selecting the same set point, yields the same concentrations again and again.

In the lab the model 3073 is ideally suited for calibration setups where low but stable aerosol concentrations are needed. Particle measurement instruments can typically be connected without additional dilution. Due to the new nozzle operating mode also high concentrations can be produced stable over long periods (VDI 3491-2).

For field use the model 3073 can be equipped with a user-supplied battery pack. This transforms the model 3073 into a completely stand-alone unit. The significantly reduced power consumption allows run times up to several hours depending on pack capacity and operating conditions. The model 3073 also features an external control interface that enables remote control which is very helpful in duct measurements.

### Features and Benefits

- Highly portable, self-contained design
- Steplessly variable output flow rate 0.3 to 4.5 L/min with improved accuracy and reliability
- Suitable for PSL, different oils, and salt solutions
- Interface supports remote control
- Internal brushless pump for long-life
- Compact and lightweight design with integrated wide-range power supply
- Optional battery for increased flexibility and field use

### Applications

- Calibration of particle sizing instruments
- General purpose reproducible test aerosol generation
- Tracer particle generation (seed particles)
- HEPA/ULPA filter integrity testing
- Cleanroom and safety cabinet validation according to ISO 14644-3
- In-use acceptance testing of low- and high-flow filters
- Smoke detector testing



## Specifications

# Portable Test Aerosol Generator

Models 3073

### Mode of Operation

Atomizer with submerged two-stream injection nozzle (Laskin mode)

### Particle Size, $d_{mode}$

0.15 to 0.3  $\mu\text{m}$  (DEHS), 0.05 to 0.2  $\mu\text{m}$  (salt, varies with concentration)

### Particle Concentration

85/cm<sup>3</sup> to >10<sup>7</sup>/cm<sup>3</sup>

### Particle Type

Liquid DEHS, PAQ (Emery 3004), Paraffin

Solid PSL, NaCl and other salt solutions

### Flow Rate

0.3 to 4.5 L/min

### Dosing Period

~44 h continuous operation

Several hours in battery operation depending on pack capacity and operating conditions

### Reservoir

100 mL standard laboratory glass bottle with DIN thread (d = 56 mm)

### Aerosol Outlet

Ø 8 mm

### Maximum Counter Pressure

40 kPa (0.4 bar)

### Dimensions (L x W x H)

30 cm x 12 cm x 19.5 cm (11.8 in. x 4.7 in. x 7.7 in.)

### Weight

3.9 kg (8.6 Lb)

### Power Supply

100-240 VAC, 50/60 Hz with IEC C13 line cable


User-supplied battery pack should provide 12 VDC and feature a 5.5 mm x 2.1 mm barrel jack for the 3073 to connect to

### Operation


The model 3073 uses a low-noise brushless swing piston pump and a submerged two-stream nozzle to produce a steady stream of aerosol. The nozzle inlet pressure can be set on the display and controls the flow and the particle production rate. As an added value of this design, the same output is achieved just by setting the nozzle pressure to the same value, even on the next day provided the solution remains the same.

As the air passes through the nozzle the Bernoulli Effect causes the liquid from which the aerosol is to be produced to be drawn in. The high shear forces within the nozzle lead to droplet formation. Since the nozzle is submerged ("Laskin mode") the droplets pass through the remaining liquid which acts as a separator for the big droplets. This defines the resulting particle size distribution in the submicrometer range.

With this design crystalline deposits are prevented and small as well as large mass flows can be generated. The fine droplet aerosol leaves the model 3073 on the top through an Ø 8 mm outlet.

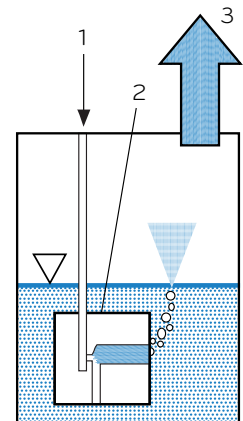

Schematic of the model 3073  
1. HEPA-filtered air flow  
2. Submerged two-stream nozzle  
3. Aerosol output

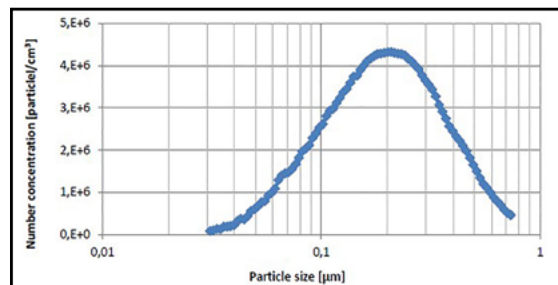

Example: Model 3073 generated DEHS aerosol measured with TSI® Scanning Mobility Particle Sizer™ (SMPS™) 3938

## To Order

### Portable Test Aerosol Generator

Specify	Description
3073	Portable Test Aerosol Generator
3062-NC/3062	Diffusion Dryer
3073-2222	External Battery Pack Connection Cable
3073-1111	Remote Control External Interface Cable

Specifications are subject to change without notice.

TSI, and the TSI logo are registered trademarks of TSI Incorporated in the United States and may be protected under other country's trademark registrations.



Knowledge Beyond Measure.

TSI Incorporated - Visit our website [www.tsi.com](http://www.tsi.com) for more information.

USA	Tel: +1 800 874 2811	India	Tel: +91 80 67877200
UK	Tel: +44 149 4 459200	China	Tel: +86 10 8219 7688
France	Tel: +33 1 41 19 21 99	Singapore	Tel: +65 6595 6388
Germany	Tel: +49 241 523030		