




Backward Compatibility:

FitPro™ Ultra Fit Test Software is Now Available for PortaCount® Fit Tester 8030 and 8038 Instruments

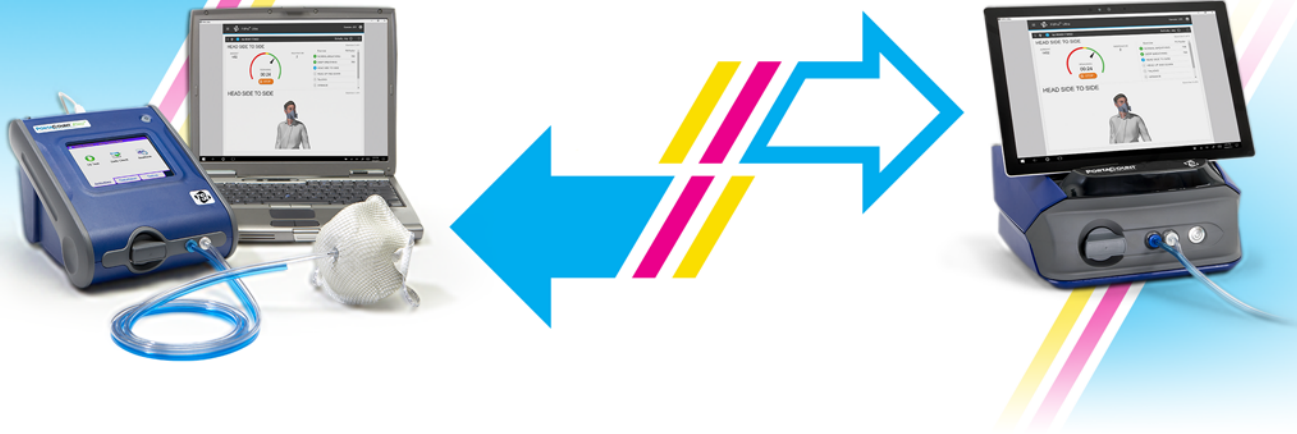

ALL OF YOUR INSTRUMENTS WORKING TOGETHER

The newest generation PortaCount® Respirator Fit Tester (models 8040 and 8048) launched with new and improved software called FitPro™ Ultra Fit Test Software (or FPU). Until now, the upgraded software was not compatible with older instruments (models 8030 or 8038).

Backward compatibility makes transitioning your program to new instruments, or working with both generations of instruments easy. Customers can now own a mix of instruments and use a single and newer software platform with connected databases.

Fit test records for older generation instruments can now reside in the newer FitPro Ultra database structure. This eliminates the need to manage and maintain two separate databases when using two generations of instruments. Backward compatibility simplifies your program, as all of your instruments use the same software.

TWO SOFTWARE VERSIONS TO SUIT YOUR NEEDS

A software license ID is required to activate the feature within the software update. Two license types are available for purchase, a basic and full version.

The basic version includes everything from FitPro Ultra, except animated fit test steps and real-time gauge features. The full version includes all software features of FitPro Ultra for your older generation instruments.

SOFTWARE COMPARISON

	Basic Version	Full Version
Core FitPro Ultra (FPU) features now available for 8030 and 8038 units		
Intelligent Fit Testing navigation, data structures and databases, reports, all other core tools	+	+
Access to future feature updates	+	+
Advanced FitPro Ultra (FPU) features now available for 8030 and 8038 units		
Fit test step animations	—	+
Real-time fit factor gauge shown during fit tests	—	+
FitCheck™ Mode with real-time gauge and graph	—	+
FitPro Ultra (FPU) features not available for 8030 or 8038 units with either license type		
Group fit testing / FitPro Ultra Client	—	—
Data import/export to/from a thumb-drive to run 8030 or 8038s in stand-alone mode	—	—



Purchase and Post-Purchase Details



NOTES ON LICENSES AND SOFTWARE

- + Licenses are limited to one computer only
 - + Once installed, you can connect to as many instruments as needed using this computer
- + To setup on another computer, additional licenses are needed
- + Licenses do not have expiration dates
- + This software will not provide enhancements that are hardware specific to newer generation instruments, such as WIFI connectivity, or lower ambient threshold levels

HOW TO PURCHASE

1. Contact your TSI Sales Rep - www.tsi.com/contact-us
2. Sales reps can explain the basics of FitPro Ultra (FPU), and the differences between the two license types
3. Select the type and number of licenses that will best serve your needs

HOW TO GET STARTED

1. After purchase, TSI will send you a license activation ID
2. Download the update to your computer (FitPro Ultra v4.8)
www.tsi.com/support/tsi-software-and-firmware-wizard
3. Follow the download wizard instructions
4. Connect your instrument to your internet-connected computer, key in your activation ID when prompted
5. Transfer your database to the new software
6. Begin using all of your instruments, working together on the same software platform and databases

If you have questions, TSI Technical Services can help you through the process. Go to **www.tsi.com/support**



UNDERSTANDING, ACCELERATED

TSI Incorporated - Visit our website **www.tsi.com** for more information.

USA Tel: +1 800 874 2811
UK Tel: +44 149 4 459200
France Tel: +33 1 41 19 21 99
Germany Tel: +49 241 523030

India Tel: +91 80 67877200
China Tel: +86 10 8219 7688
Singapore Tel: +65 6595 6388