

CONDENSATION PARTICLE COUNTER MODEL 3750

THE MOST REFERENCED CPCs
WITH OVER 40 YEARS
OF EXPERIENCE

The 4th generation of TSI Butanol-CPCs (Condensation Particle Counters) makes the best CPCs even better: providing easier use for everyone with an intuitive direct-touch user-interface, new software and secure internal data storage for months of data. Users can rely on the proven compact and rugged instrument with new wired or optional wireless communication capabilities accommodating the modern needs of researchers, enabling full remote access to instruments located in monitoring sites around the world.



Features and Benefits

- + Extended concentration range up to 100,000 particles/cm³
- + Extended size range down to 7 nm (D50) -
Improved counting efficiency for SMPS¹
- + Data rate up to 50 Hz -
Capture highly dynamic processes
- + Water-removal system -
Compatible to high-humidity environments
- + Improved Diagnostics with Pulse Height Monitoring
- + Upgradeable with CEN/TS 16976 calibration for continuous ambient monitoring
- + Upgradeable to 1-nm-CPC with compatible Nano Enhancer


The Standard CPC model 3750 fits seamlessly into the family of four Butanol-based CPCs including the High Concentration, the Ultrafine Particle and the 1-nm model, and is fully integrated into the SMPS¹ family series 3938. With its high counting flow, concentration and size range it serves all applications from filter testing, basic aerosol research, combustion, health and atmospheric studies.

All TSI CPCs are equipped with a removable saturator wick for easy transport and maintenance to allow longer operation time by simple user maintenance.

¹ SMPS stands for Scanning Mobility Sizer, another reference system by TSI to measure nanoparticle size distributions.



UNDERSTANDING, ACCELERATED

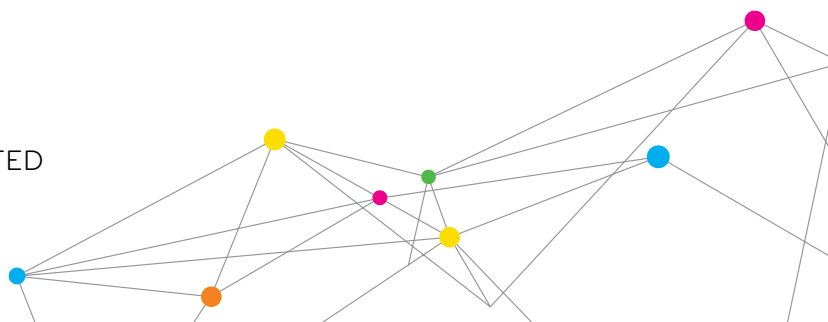

SPECIFICATIONS

CONDENSATION PARTICLE COUNTER MODEL 3750

Particle Size Range

7 nm min. detectable particle size (D50), verified with monodisperse sucrose particles
>3 µm max. detectable particle size

Particle Concentration Range

Up to 100,000 (1x10⁵) particles/cm³
Single particle counting with continuous live-time coincidence correction

Particle Concentration Accuracy

±5% at <100,000 particles/cm³

False Background Counts

<0.001 particles/cm³ based on 12 hour average

Response Time

Response time described in percentage to concentration step change
<1 second for 90% to 10% (T10-90, T90-10)
~2 seconds for 0 to 95% (T95)

Flow System

1.0 ± 0.05 L/min inlet and counting flow (volumetric)
Requires external pump

Liquid System

Butanol (n-Butyl alcohol, not included) used as working fluid
Use of the internal water removal pump to remove condensate is recommended for use in humid environments

Communication Interfaces

Embedded touch-display
Pulse output: BNC connector, TTL level pulse, nominally 350 nanoseconds wide
USB type C to connect CPC directly to computer operating control software Aerosol Instrument Manager (included).
Ethernet port (8-wire RJ-45 jack, 10/100 BASE-T, TCP/IP) for remote connection. Configurable for automated (DHCP) or manual network settings.
USB type B for external memory drives.
A Wi-Fi adapter can be used to enable wireless network integration.

Ambient Operating Conditions

Temperature 10 to 35°C (50 to 95°F)
Humidity 0 to 90% RH, noncondensing
Pressure 75 to 105 kPa (0.75 to 1.05 atm)

Accessories

Required

Electrical: 100 to 240 VAC, 50/60 Hz, 200 W maximum.
Auto recovery from power failure built in
Vacuum: 60 kPa (18 in Hg) minimum gauge (below atmospheric pressure)

Included

Fill and drain bottles
Aerosol Instrument Manager for Count products license
USB C to A cable for connecting to a computer

Data Storage

Internal memory lasts for approx. 1 year of data at 50 Hz data rate.

Dimensions (H x W x D)

27.5 cm x 18.3 cm x 29.9 cm (10.83 in. x 7.21 in. x 11.76 in.), not including fill bottle and bracket

Weight

~6.6 kg (~14.6 lbs)

Catalog Items

Specify	Description
3750	Condensation Particle Counter
3032	Vacuum pump 110V
3032-1	Vacuum pump 230 V/50 Hz
3750-	MKIT Maintenance Kit for CPC
3750-WKIT	Wick replacement kit for CPC
3772200	Environmental sampling system

Specifications are subject to change without notice.

TSI and the TSI logo are registered trademarks of TSI Incorporated.



UNDERSTANDING, ACCELERATED

TSI Incorporated - Visit our website www.tsi.com for more information.

USA Tel: +1 800 874 2811 **India** Tel: +91 80 67877200
UK Tel: +44 149 4 459200 **China** Tel: +86 10 8219 7688
France Tel: +33 1 41 19 21 99 **Singapore** Tel: +65 6595 6388
Germany Tel: +49 241 523030