

Dilution System

Model 3332




Application Note 3332-001 (US)

Dilution System for 1 Liter Per Minute Flows

The main application of the Model 3332 Dilution System is in the dilution of test aerosols for the Model 3330 Optical Particle Sizer. Other applications having a 1 liter per minute flow rate can also use the Model 3332 Particle Diluter.

The diluter is available in two models with different dilution factors.

	Volume Flow	Dilution Factor
3332-100	1.0 L/min	100

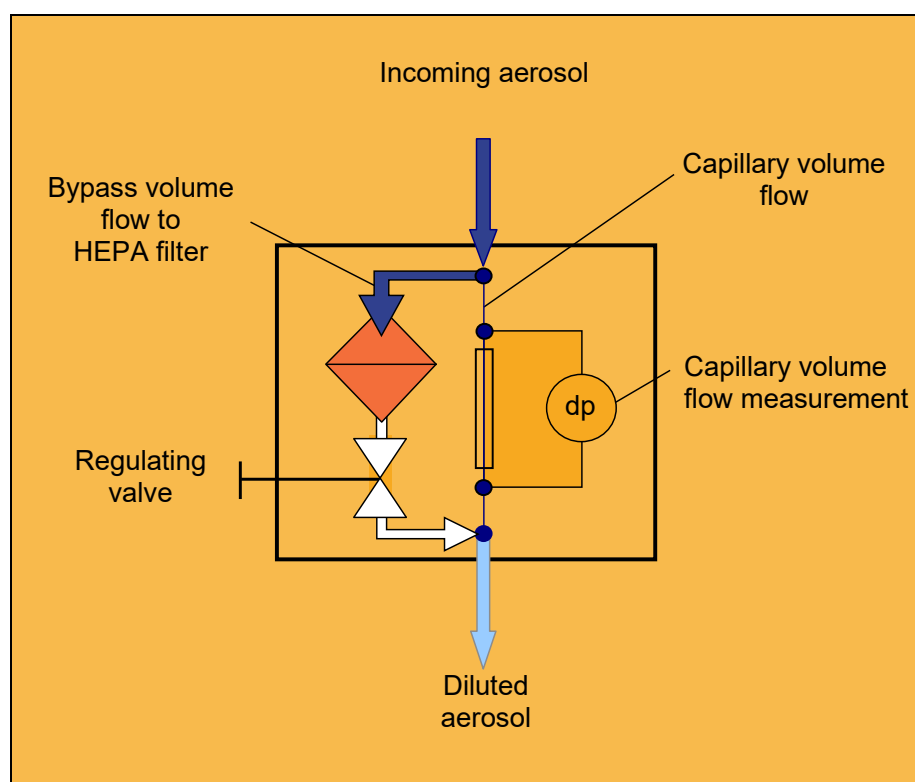

Schematic of the Model 3332 Dilution System



Model 3332 Dilution System

The indicator lights on the diluter indicate if the pressure drop across the capillary is what it should be (based on its calibration). If one of the red arrows is illuminated, an adjustment of the knob will restore the proper balance between the flows.

Power Supply	9 ... 15 VDC (AC adapter, battery operation)
Dimensions	160 mm x 300 mm x 110 mm (6.3 in. x 11.8 in. x 4.3 in)

TSI and the TSI logo are registered trademarks of TSI Incorporated in the United States and may be protected under other country's trademark registrations.



Knowledge Beyond Measure.

TSI Incorporated – Visit our website www.tsi.com for more information.

USA Tel: +1 800 680 1220
UK Tel: +44 149 4 459200
France Tel: +33 1 41 19 21 99
Germany Tel: +49 241 523030

India Tel: +91 80 67877200
China Tel: +86 10 8219 7688
Singapore Tel: +65 6595 6388