


Service Options



Laser Aerosol Spectrometer 3340A

TSI® service, performed on an annual basis, helps to keep your LAS 3340A in proper working order. TSI® offers two levels of service: Clean and Calibrate, and Repair, Clean and Calibrate. Read below to learn what is included in each service.

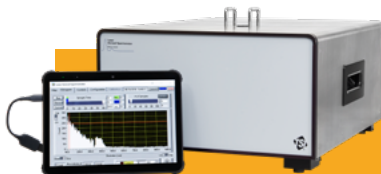

Clean and Calibrate Service includes:

Maintain

- Confirm voltages at circuit boards. Set reference voltage to ensure proper sizing and counting
- Verify fan performance and leak-free flow path
- Clean laser window
- Align optics while running reference particles (100 nm PSL)

Replace

- Replace filters, tubing, and o-rings as needed

- Calibrate particle sizing from 0.1 – 7 microns using NIST-traceable PSL

Zero Check

Test

- Conduct 24 hours of zero check testing to ensure concentration accuracy, and meet or exceed false count specification of < 1 particle / 5 minutes

Counting Efficiency

- Confirm counting efficiency of instrument relative to a reference single-particle-counting instrument (CPC) using monodisperse aerosol

Does NOT Include


Broken Parts

- Internal pump, circuit board, etc.
(In such cases, consider 'Repair, Clean and Calibrate' service)

Physical Damage

- Scratches on case

Service Options Laser Aerosol Spectrometer 3340A

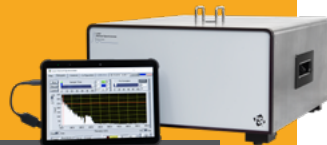

Repair, Clean and Calibrate includes:

Repair • Repair or replace broken parts (e.g. internal pump, circuit board)

Clean and Calibrate • All actions in a standard 'Clean and Calibrate' service are then performed

Does NOT Include • Physical damage (e.g. scratches on case)

What Makes TSI® Service Valuable?

As a result of TSI® service, you are returned an optimally performing and updated instrument. The service is performed by TSI's experienced technicians for a fixed price. This simple service process allows you to stay focused on the work that is important to you. TSI's calibration of the 3340A LAS is traceable to National Standards and Metrology Institutes (NIST, NPL UK) and other accepted values of physical constants. TSI® is registered to ISO-9001:2015.

These TSI services come with a 90-day warranty on labor performed, or parts used, during the service.

Request Service

Specify	Description
CL-3340A	Clean and Calibrate
RP-3340A	Repair, Clean and Calibrate

Call, or visit
tsi.com/requestservice



Knowledge Beyond Measure.

TSI Incorporated - Visit our website www.tsi.com for more information.

USA Tel: +1 800 874 2811
UK Tel: +44 149 4 459200
France Tel: +33 1 41 19 21 99
Germany Tel: +49 241 523030

India Tel: +91 80 67877200
China Tel: +86 10 8219 7688
Singapore Tel: +65 6595 6388